

PAWLENX · 2026-05-17

UI/UX Flow Diagrams


Mobile app product flows · annotated for design hand-off

14 diagrams · each followed by a description page that names every node.

Contents

Each diagram is on its own large page; a description page follows it where node labels need explanation.

1.	Top-level user journey	8.	Companion AI Chat (F-CHAT)
2.	Information architecture (sitemap)	9.	Preventive care & reminders (F-CARE)
3.	Diagnosis & triage flow (flagship — F-DIAG)	10.	Vet handoff (F-VET)
4.	Onboarding & first pet profile	11.	Treatment & recovery (F-RX)
5.	Home / Daily Companion	12.	Household & multi-pet (F-HH)
6.	Timeline & trends (F-TL)	13.	Cross-cutting AI trust pattern
7.	Nutrition engine (F-NUT)	14.	State matrix per screen (required deliverables)

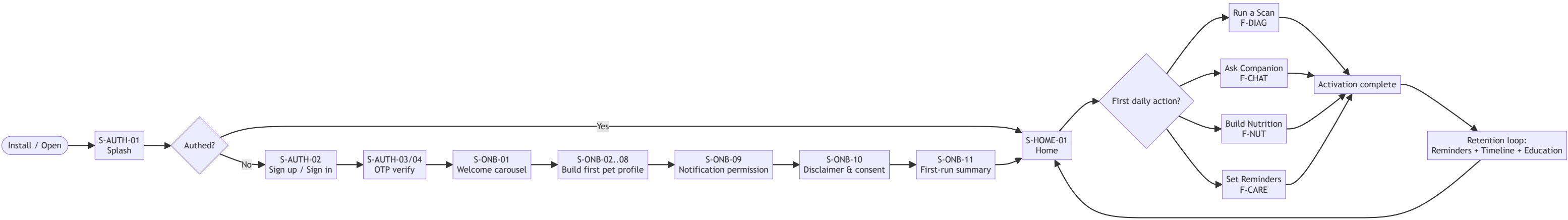

Top-level user journey

Acquisition → activation → retention loop. Shows how a new user moves from install through onboarding into their first daily action and back into the retention loop.

ENTRY & AUTH

Install / Open

User launches the app for the first time, or returns to it.

S-AUTH-01 Splash

Branded splash while the app boots, restores session and decides where to route.

Authenticated?

Decision: is there a valid session? Yes → straight to Home. No → enter sign-up/sign-in.

S-AUTH-02 Sign up / Sign in

Email + Apple/Google. Handles both new and returning users.

S-AUTH-03/04 OTP verify

One-time code verification screens for email or phone.

ONBOARDING

S-ONB-01 Welcome carousel

3–4 swipeable cards explaining the core promise (diagnose, plan, track).

S-ONB-02..08 Build first pet profile

Sequenced steps to capture species, name, breed, weight, lifestyle, history and a photo.

S-ONB-09 Notification permission

OS push prompt framed with a "why we ask" rationale.

S-ONB-10 Disclaimer & consent

Plain-language statement that PawLenx is decision support, not a vet. Required to proceed.

S-ONB-11 First-run summary

Suggests the user’s very first action (scan, ask, plan, remind).

HOME & ACTIVATION

S-HOME-01 Home

Default landing screen — Today card, quick actions, feed, FAB.

First daily action?

Branch tracked by analytics: the first meaningful action determines activation type.

Run a Scan · Ask Companion · Build Nutrition · Set Reminders


Four primary activation paths, each linked to a feature ID (F-DIAG, F-CHAT, F-NUT, F-CARE).

Activation complete

User has done one meaningful action; flagged in analytics as activated.

Retention loop

Reminders, timeline updates and education cards pull the user back to Home over time.

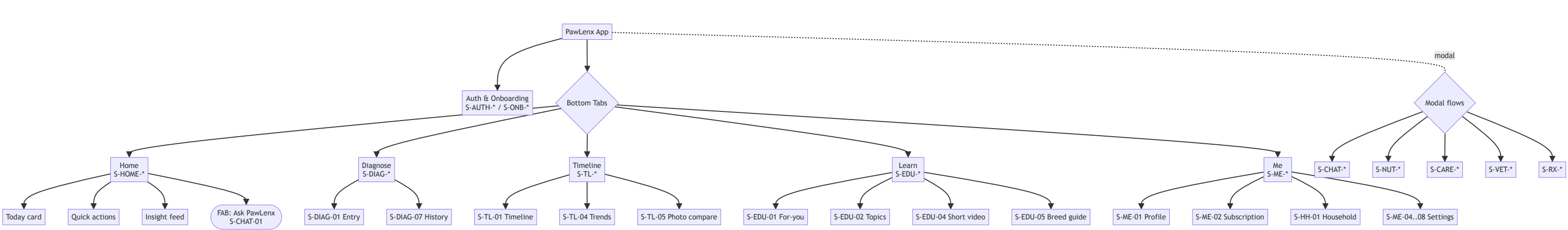

Information architecture (sitemap)

Top-level navigation. Five bottom tabs are the persistent surface; modal flows are summoned from tabs or from the floating "Ask PawLenx" button.

BOTTOM TABS

Home (S-HOME-*)

Personalised daily companion; first thing users see after auth.

Diagnose (S-DIAG-*)

Entry point to the scan/triage feature plus scan history.

Timeline (S-TL-*)

Chronological pet health record + trend charts + photo comparisons.

Learn (S-EDU-*)

For-you feed, topic library, short video and breed guides.

Me (S-ME-*)

Profile, subscription, household management, settings.

MODAL FLOWS

S-CHAT-*

Companion AI chat, opened from the global FAB or from any AI result.

S-NUT-*

Nutrition planning, food picker, allergen check.

S-CARE-*


Preventive care plan and reminders.

S-VET-*

Vet handoff: auto-summary, edit, share or export PDF.

S-RX-*

Active treatments, dose reminders, recovery check-ins.

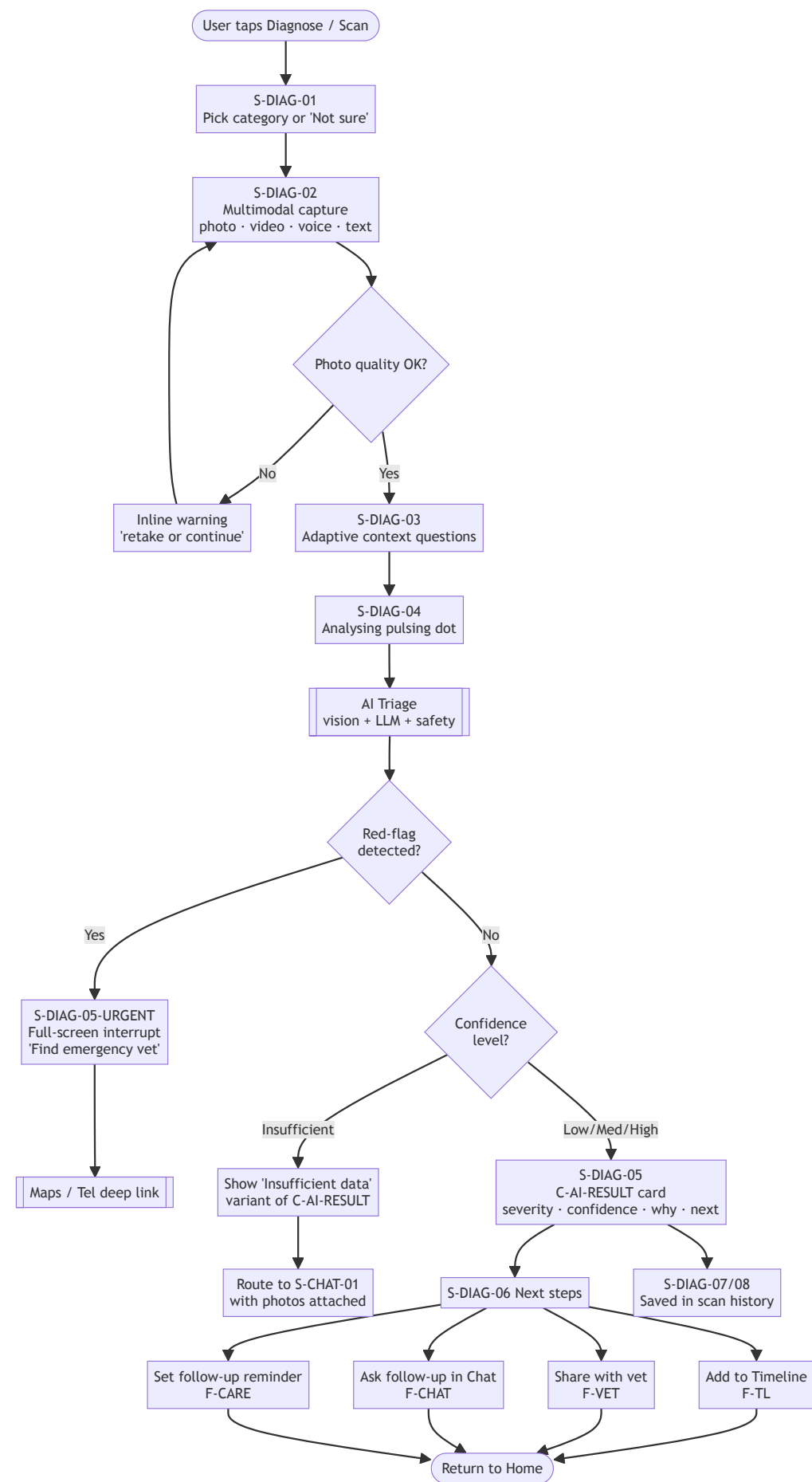

Diagnosis & triage flow (flagship — F-DIAG)

The most important flow in the product. Mandatory red-flag interrupt; confidence-aware results; every result is also stored as a timeline event.

CAPTURE

S-DIAG-01 Pick category or "Not sure"

User selects a body system (skin, eyes, GI...) or "Not sure" to let AI infer.

S-DIAG-02 Multimodal capture

Camera, video, voice and text input on one screen with on-device quality guidance.

Photo quality OK?

On-device check (blur, lighting, framing). Failure shows inline retake hint.

S-DIAG-03 Adaptive context questions

Dynamic 3–6 question set based on category; answers feed the AI.

AI RUN & SAFETY

S-DIAG-04 Analysing

Pulsing-dot waiting state with progress hint; cannot navigate away.

AI Triage

Vision classifiers + grounded LLM + safety pipeline (see Engineering #5).

Red-flag detected?

Hard gate. If yes, the normal result UI is replaced by an emergency interrupt.

S-DIAG-05-URGENT

Full-screen, non-dismissable interrupt: "Find emergency vet" + tel/maps deep links. No upsell.

Confidence level?

Tri-state: Insufficient (route to chat), Low/Medium/High (show result with confidence band).

Insufficient data

Variant of the result card that admits uncertainty and offers to continue in chat.

RESULT & NEXT STEPS

S-DIAG-05 C-AI-RESULT card


Severity, confidence dots, "why we think this", next steps, citations, disclaimer.

S-DIAG-06 Next steps

Four primary CTAs: timeline, vet, chat follow-up, reminder.

S-DIAG-07/08 Scan history

Every scan is persisted and listed under the Diagnose tab.

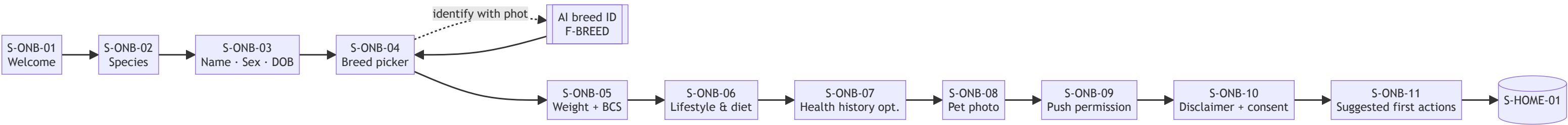

Onboarding & first pet profile

A guided, skippable wizard that captures enough information for the AI to be useful from day one. Each step has its own screen ID.

STEPS

S-ONB-01 Welcome

Brand carousel; sets expectations.

S-ONB-02 Species

Dog / Cat (more in V2).

S-ONB-03 Name · Sex · DOB

Basic identity. DOB drives life-stage logic.

S-ONB-04 Breed picker

Searchable list. Optional photo → AI breed ID (F-BREED).

S-ONB-05 Weight + BCS

Weight in kg/lb plus a 5-image Body Condition Score chooser.

S-ONB-06 Lifestyle & diet

Indoor/outdoor, activity, current food brand/type.

S-ONB-07 Health history (optional)

Allergies, chronic conditions, vaccinations.

S-ONB-08 Pet photo

Cover photo for the profile and timeline.

S-ONB-09 Push permission


OS prompt for reminders.

S-ONB-10 Disclaimer + consent

Required legal screen.

S-ONB-11 Suggested first actions

Personalised "do this next" prompts.

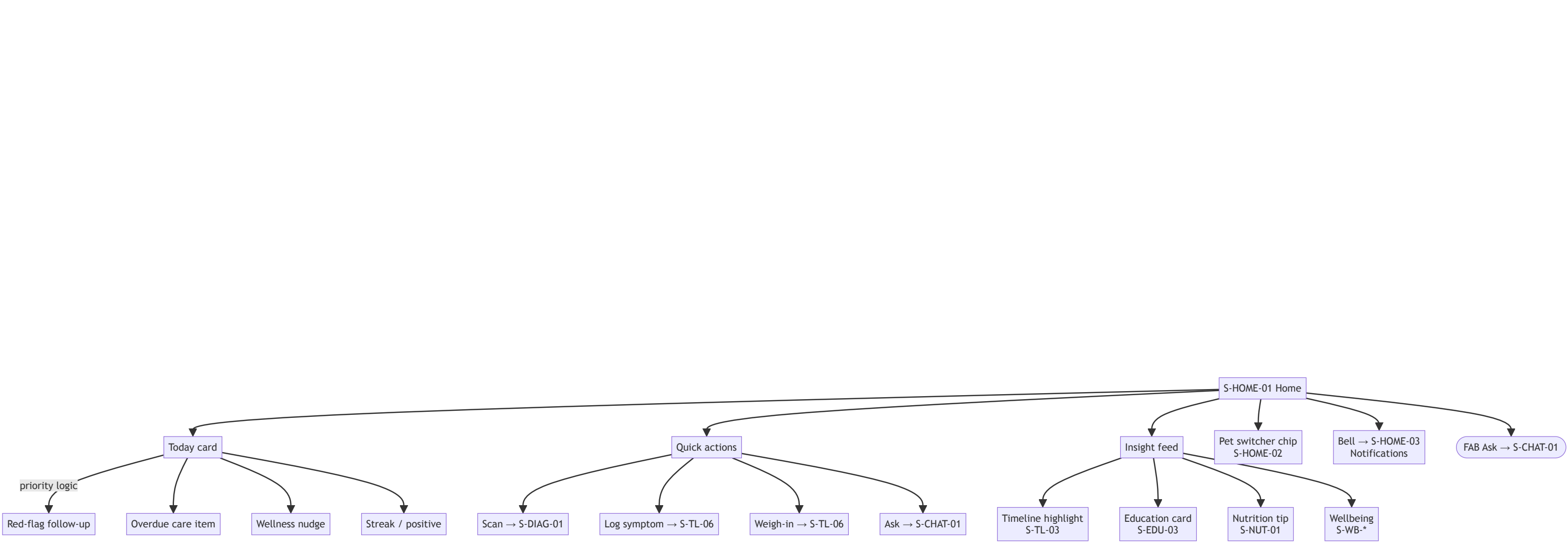

Home / Daily Companion

A single screen that picks the most important thing for today and offers the four most-used actions plus a personalised feed.

TODAY CARD (ONE PRIORITY CHOSEN BY RANKER)


- Red-flag follow-up**
Highest priority — yesterday’s urgent result needs an outcome.
- Overdue care item**
Missed vaccination, medication or recurring task.
- Wellness nudge**
Gentle, evidence-based suggestion (e.g. weigh-in).
- Streak / positive**
Celebratory message when nothing else is urgent.

QUICK ACTIONS

- Scan → S-DIAG-01**
One tap into the diagnose flow.
- Log symptom → S-TL-06**
Quick free-text or category log.
- Weigh-in → S-TL-06**
Recording weight feeds trend charts.
- Ask → S-CHAT-01**
Free-form question to the AI companion.

FEED & GLOBAL

- Insight feed**
Mix of timeline highlights, education, nutrition tips, wellbeing prompts.
- Pet switcher chip (S-HOME-02)**
Households can have multiple pets; this switches active pet context.
- Bell → S-HOME-03 Notifications**
In-app notification list.
- FAB Ask → S-CHAT-01**
Persistent floating "Ask PawLenx" button.

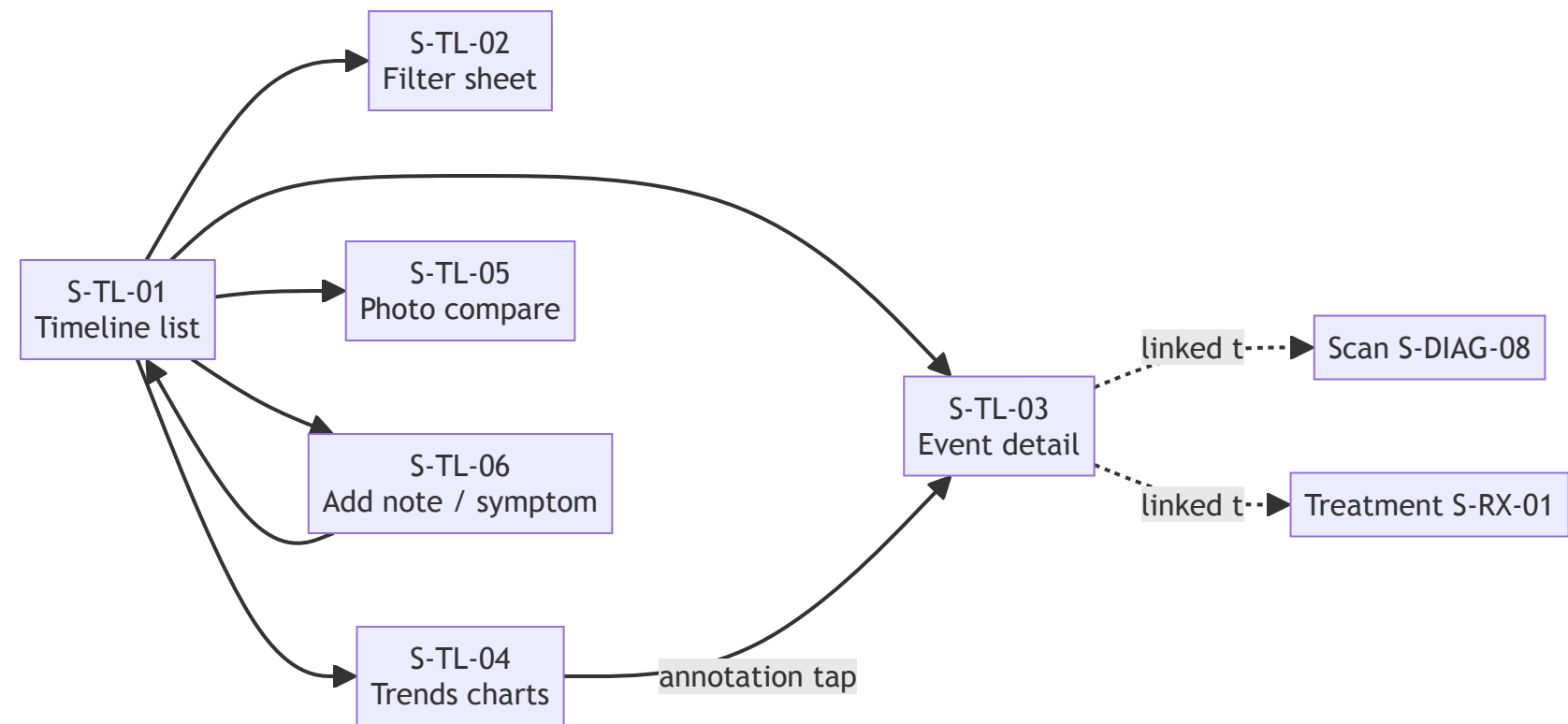

Timeline & trends (F-TL)

The pet’s health record. Every scan, treatment, note, weight and care event lives here. Designed to be the screen the user shows the vet.

SCREENS

S-TL-01 Timeline list

Reverse-chronological feed of events grouped by month.

S-TL-02 Filter sheet

Filter by type, body system, severity, date range.

S-TL-03 Event detail

Single event view with linked scan, treatment or note.

S-TL-04 Trends charts


Weight, BCS, symptom frequency over time. Tap a point to jump to that day.

S-TL-05 Photo compare

Side-by-side photos of the same body part across time.

S-TL-06 Add note / symptom

Manual entry sheet.

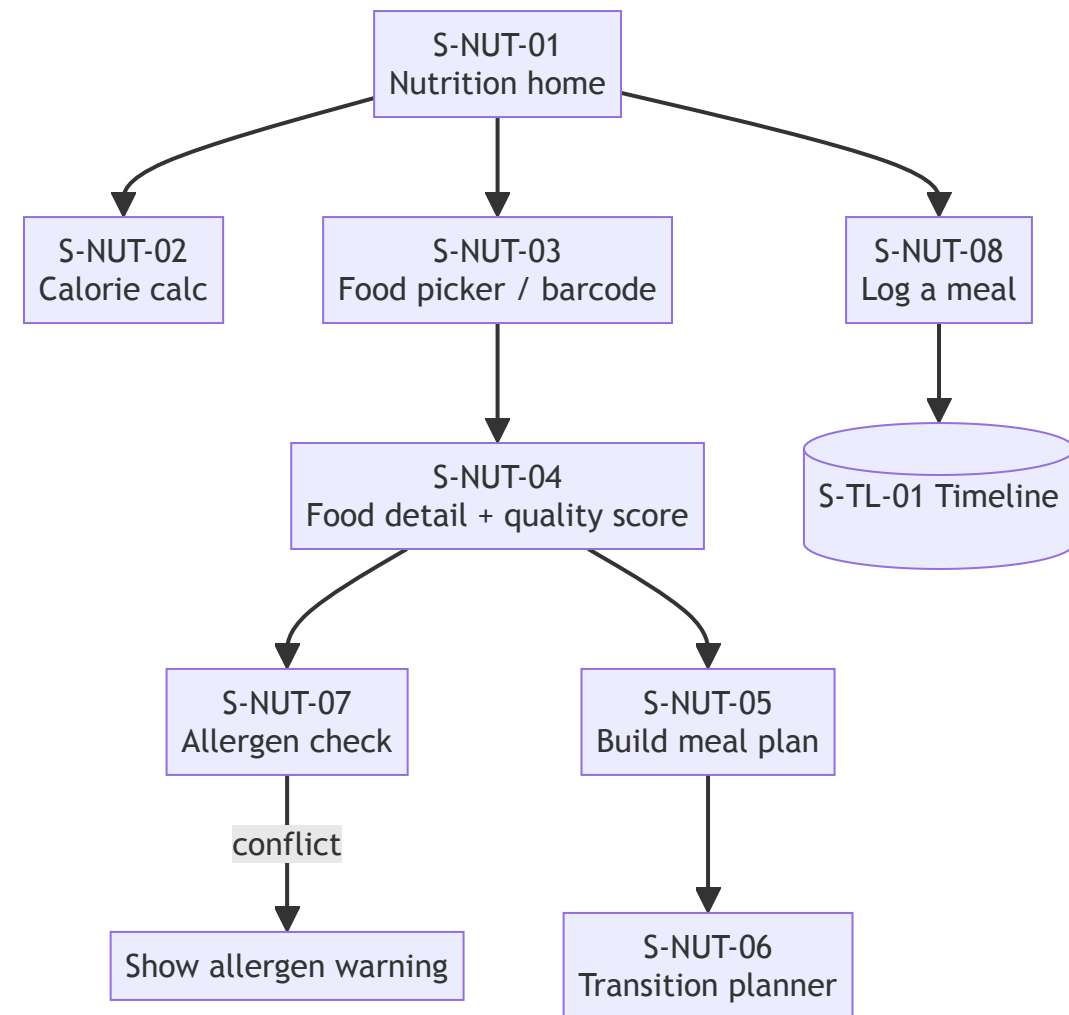

Nutrition engine (F-NUT)

Personalised feeding guidance with allergen safety and food database lookup. Logging a meal feeds the timeline.

SCREENS

S-NUT-01 Nutrition home

Plan overview, today’s meals, quick log.

S-NUT-02 Calorie calc

Daily kcal target from weight, BCS, activity, life-stage.

S-NUT-03 Food picker / barcode

Search or scan a product to add it to the plan.

S-NUT-04 Food detail + quality score

Ingredients, AAFCO/FEDIAF compliance, our quality score.

S-NUT-05 Build meal plan

Schedule meals across the day.

S-NUT-06 Transition planner


Stepped switchover schedule between two foods.

S-NUT-07 Allergen check

Cross-references the pet’s known allergies against the food.

S-NUT-08 Log a meal

Quick log → appears as a timeline event.

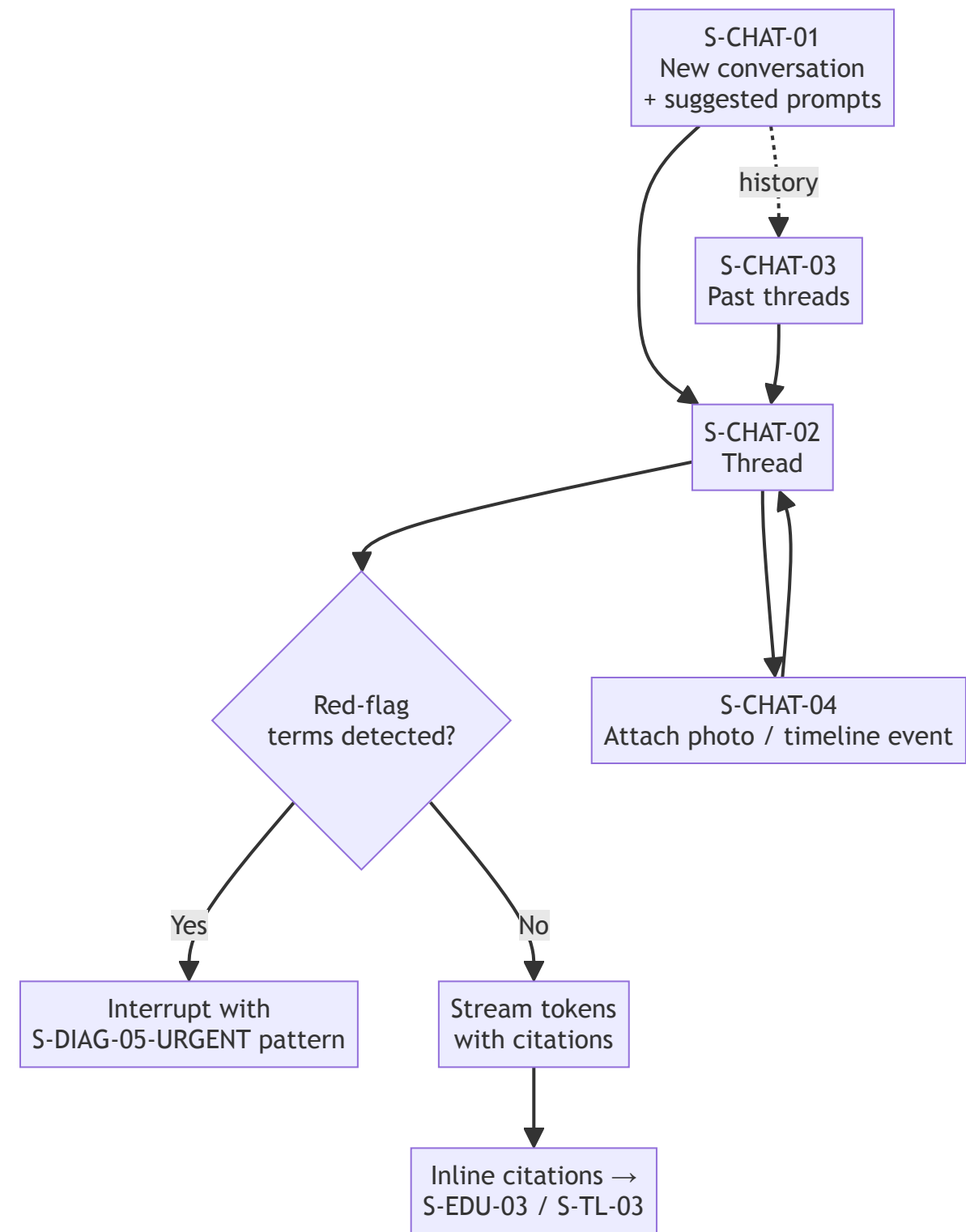

Companion AI Chat (F-CHAT)

Free-form conversation grounded in the pet’s history and the curated knowledge base. The same red-flag interrupt as the diagnose flow applies here.

SCREENS

S-CHAT-01 New conversation

Blank state with suggested prompts and recent threads.

S-CHAT-02 Thread

Streaming tokens, inline citations, attach actions.

S-CHAT-03 Past threads

Browse/search history.

S-CHAT-04 Attach photo / timeline event

Send media or pull a past event into the conversation.


BEHAVIOUR

Red-flag terms detected?

Server-side check; if hit, the chat is interrupted with the S-DIAG-05-URGENT pattern.

Stream tokens with citations

Tokens arrive over SSE; citations link to S-EDU-03 articles or S-TL-03 events.

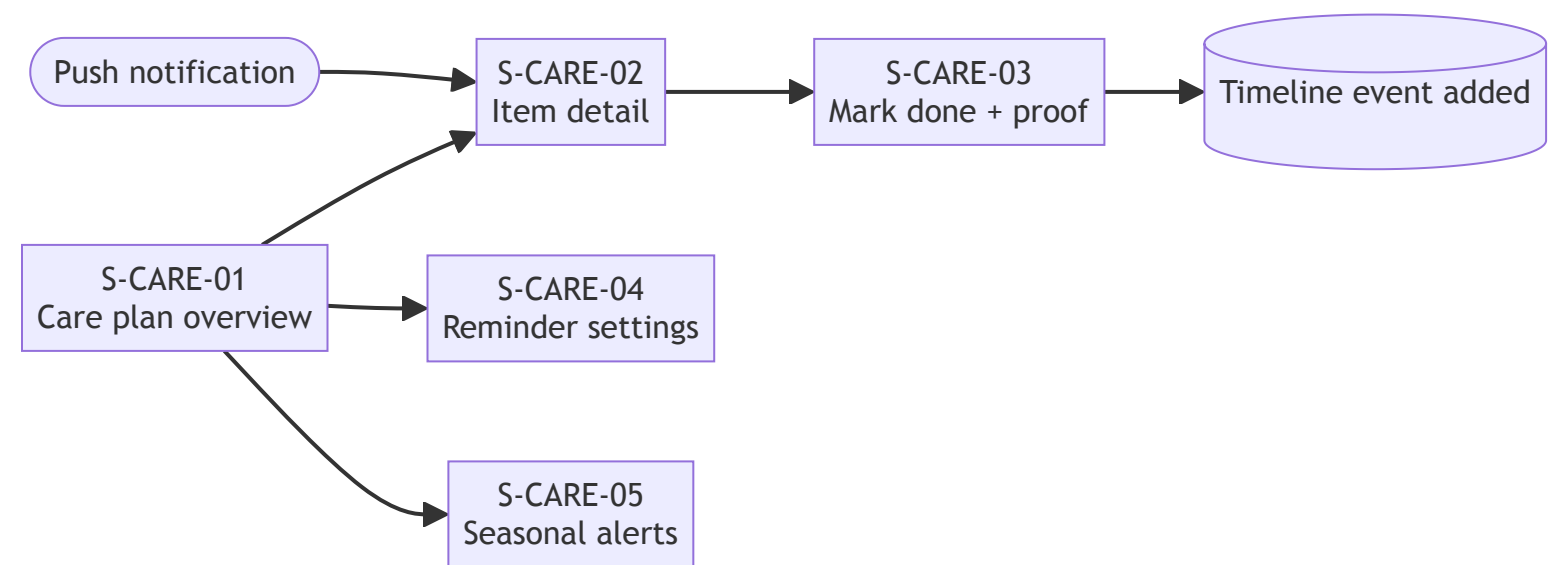

Preventive care & reminders (F-CARE)

Care plan management — vaccinations, parasite control, dental, grooming, seasonal alerts.

SCREENS

S-CARE-01 Care plan overview

List of upcoming and overdue items per pet.

S-CARE-02 Item detail

Background, recommended interval, history.

S-CARE-03 Mark done + proof

Tick off with optional photo/receipt; creates a timeline event.

S-CARE-04 Reminder settings


Per-item lead time, time of day, quiet hours.

S-CARE-05 Seasonal alerts

Region-aware (ticks, heat, holidays).

Push notification

Deep-links straight into the relevant care item.

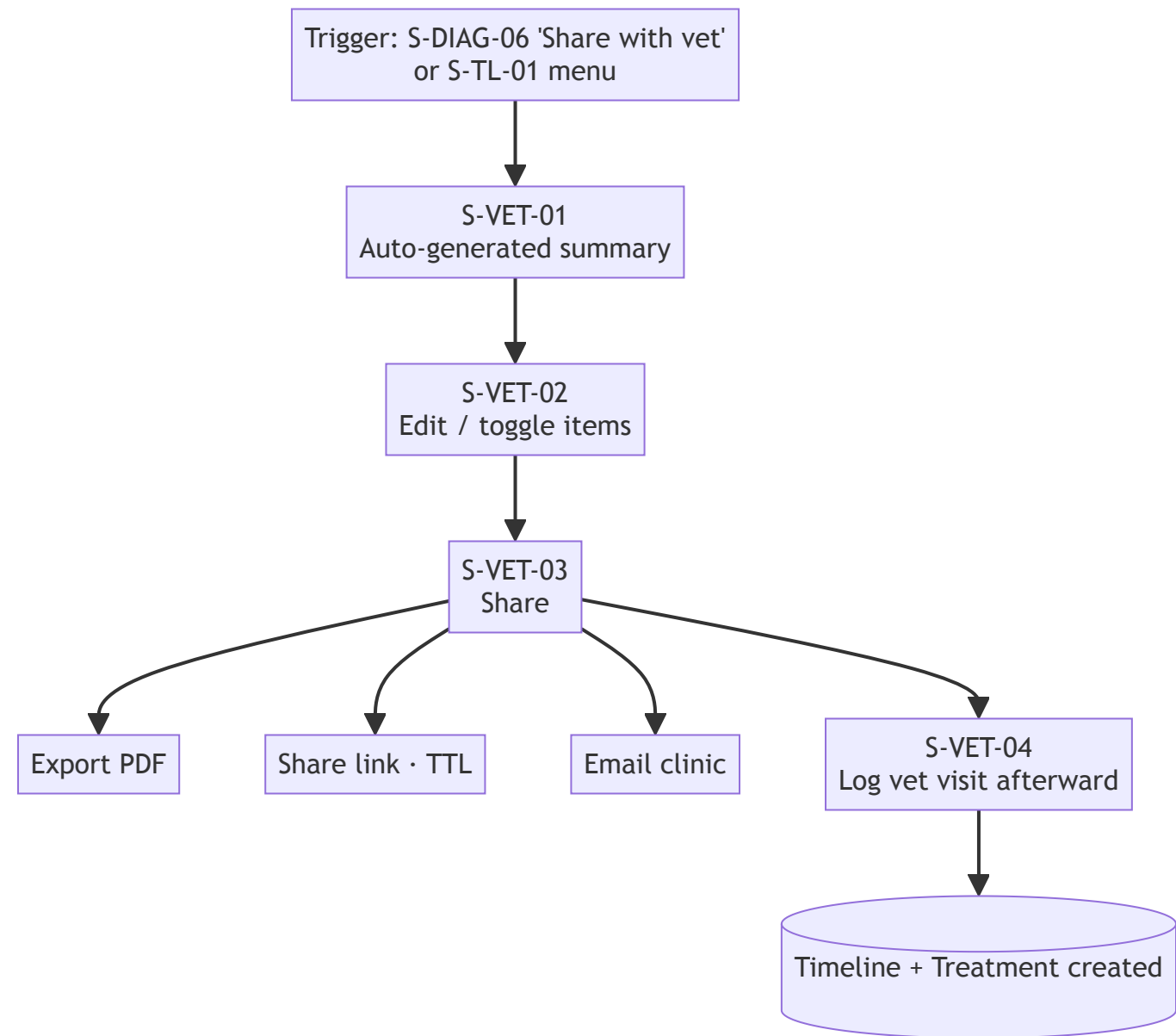

Vet handoff (F-VET)

Turn the in-app data into a vet-friendly summary in seconds — exportable as a PDF or a time-limited share link.

SCREENS

S-VET-01 Auto-generated summary

Server compiles recent symptoms, scans, weights, meds into a single document.

S-VET-02 Edit / toggle items


User can hide individual items before sharing.

S-VET-03 Share

Three sub-actions: export PDF, generate a TTL link, or email the clinic.

S-VET-04 Log vet visit afterward

Capture diagnosis, prescriptions and follow-ups; creates timeline + treatment entries.

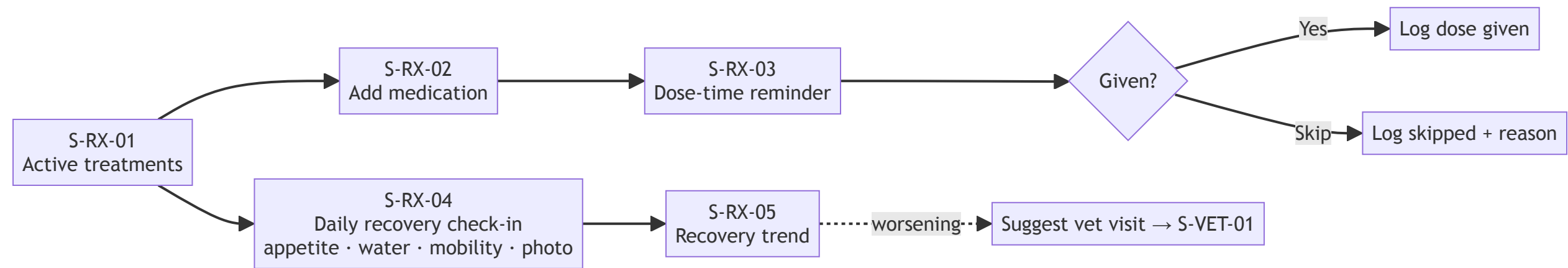

Treatment & recovery (F-RX)

Medication scheduling and daily recovery check-ins after illness or surgery.

SCREENS

S-RX-01 Active treatments

List of current medications and recovery plans.

S-RX-02 Add medication

Name, dose, frequency, duration; can be pre-filled from S-VET-04.

S-RX-03 Dose-time reminder


Push at the right time; user marks "given" or "skipped" with reason.

S-RX-04 Daily recovery check-in

Quick sliders: appetite, water, mobility + a photo.

S-RX-05 Recovery trend

Chart of the check-in metrics; worsening trend suggests a vet revisit.

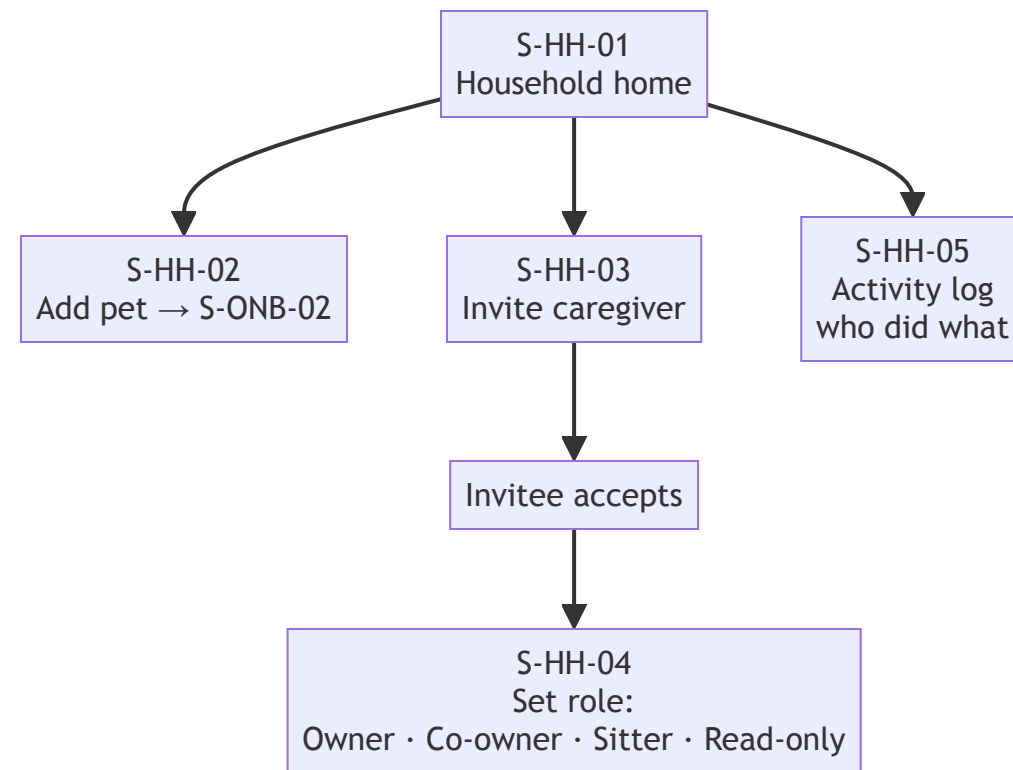

Household & multi-pet (F-HH)

One account can manage multiple pets and invite other caregivers (partner, family, sitter).

SCREENS

S-HH-01 Household home

List of pets + members in the household.

S-HH-02 Add pet

Re-uses the onboarding flow from S-ONB-02 onwards.

S-HH-03 Invite caregiver


Email/SMS invite with a pending state.

S-HH-04 Set role

Owner · Co-owner · Sitter · Read-only — controls what actions are allowed.

S-HH-05 Activity log

Audit trail of who logged what (transparency between caregivers).

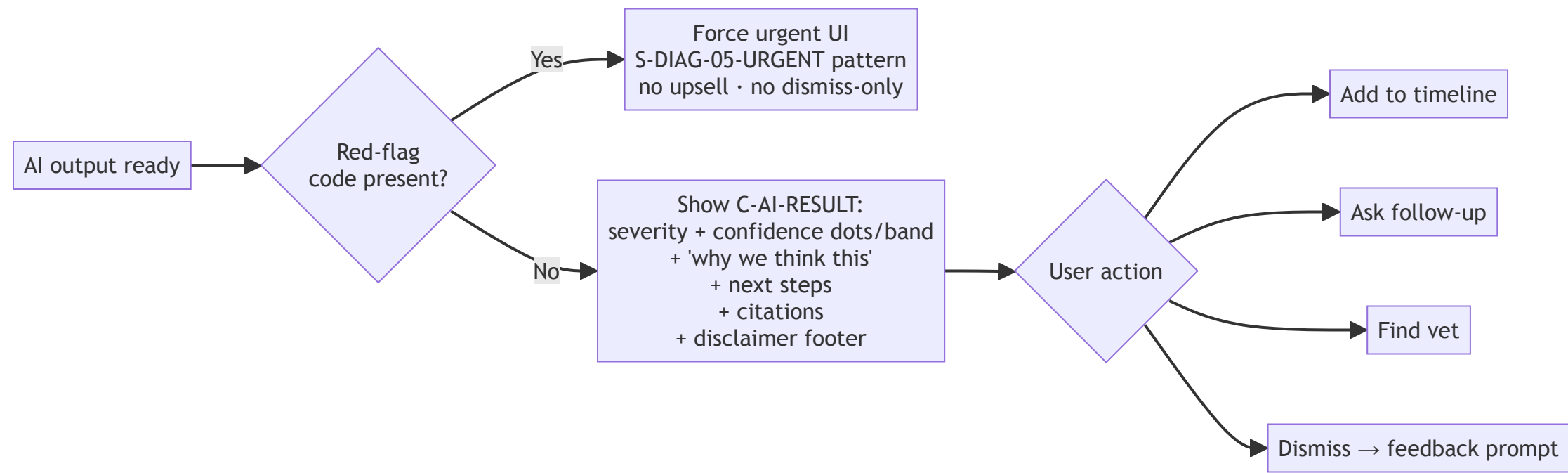

Cross-cutting AI trust pattern

Every AI output anywhere in the product passes through this UI pattern. Not optional. Diagnose, chat, nutrition, breed ID — all use it.

STATES

Red-flag code present?

Server sets one or more codes (e.g. RF-BLOAT, RF-SEIZURE). UI must surface them differently.

Force urgent UI (S-DIAG-05-URGENT)

No upsell, no dismiss-only. Phone + maps shortcuts.

C-AI-RESULT card

The standard component: severity, confidence dots/band, "why we think this", next steps, citations, disclaimer footer.

User action

Add to timeline · ask follow-up · find vet · dismiss → feedback prompt on dismiss.

